


# Blood Glucose Test

## Medi-Test Glycaemie C



- Visual blood glucose determination
- Reliable results without instrument
- No calibration necessary

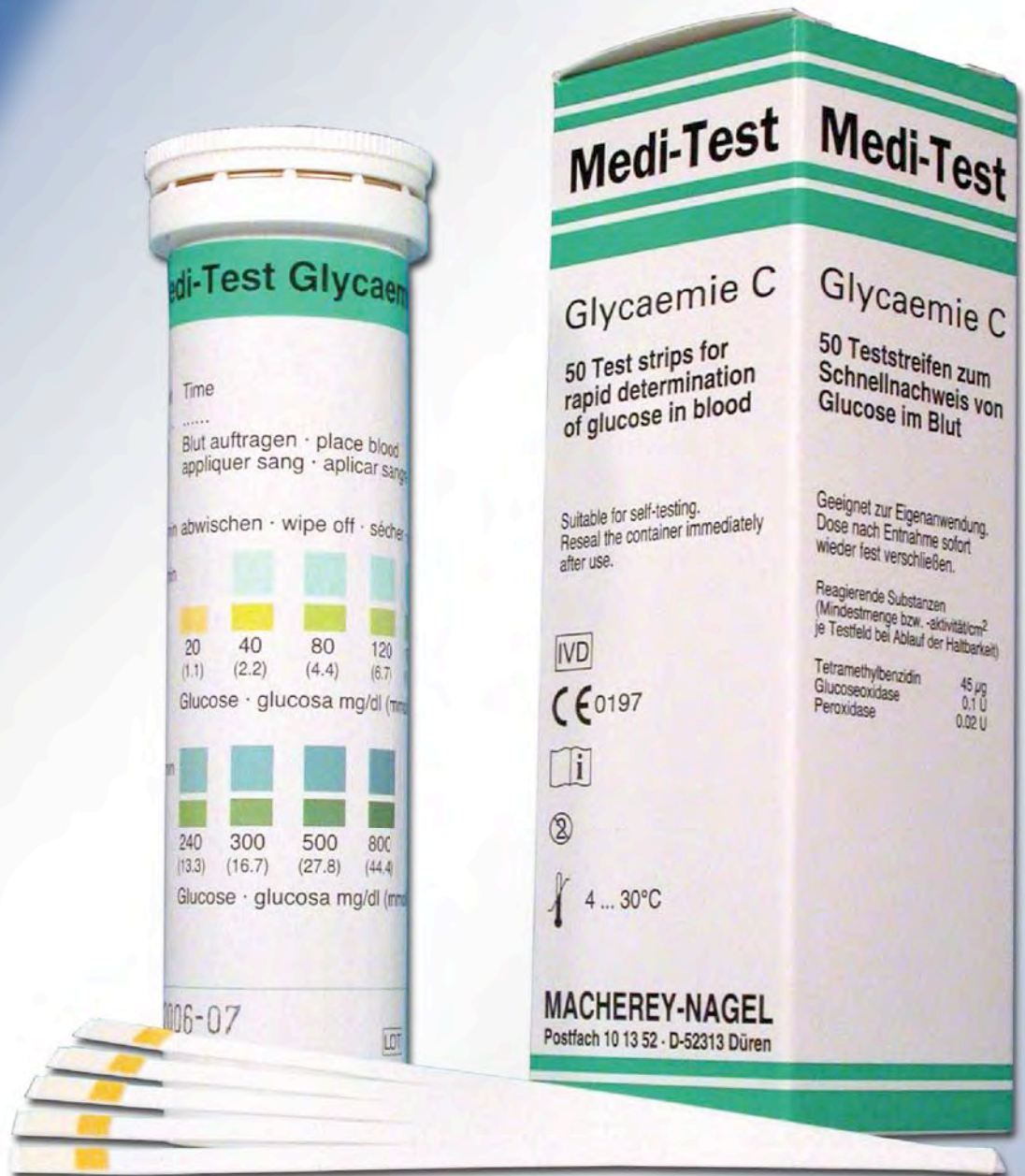

## **Always ready for use**

Strips can easily be stored and transported. Just open the container, remove one strip and in short time you get a first impression of the patients blood glucose level. Instrumental evaluation of the strip is not needed.

## **Easy measurement**

Use of Glycaemie C is easy:

- Apply one drop of blood to the test field.
- Wipe off excess blood after 1 minute.
- Read the result after another 2 - 3 minutes.

## **Best reliable results**

Two indicator zones combined with a clear colour chart ensure best reliable results wherever you are. Instruments may give false results when not correctly calibrated. The visual test strips Glycaemie C are always correct.

## **CE-marked**

The IVD-CE according to the directive 98/79/EG ensures both optimal quality and optimal usability by health care professionals as well as untrained persons.



## **Ordering information:**

Glycaemie C, blood glucose test strips  
for visual evaluation;  
gradation: 20-40-80-120-180-240-300-500-800 mg/dl;  
CE-marked, suitable for self-testing;  
box of 50 strips

**Cat. No. 928 07**

MACHEREY-NAGEL – more than 25 years experience in medical test strips

**MACHEREY-NAGEL**

MACHEREY-NAGEL GmbH & Co. KG · Postfach 10 13 52 · D-52313 Düren

Germany and worldwide:

Tel.: +49 (0) 24 21 96 90

Fax: +49 (0) 24 21 96 91 99

e-mail: sales-de@mn-net.com

Switzerland:

Tel.: +41 (0) 62 388 55 00

Fax: +41 (0) 62 388 55 05

e-mail: sales-ch@mn-net.com

France:

Tel.: +33 (0) 3 88 68 22 68

Fax: +33 (0) 3 88 51 76 88

e-mail: sales-fr@mn-net.com



[www.mn-net.com](http://www.mn-net.com)